

Roosa Kuusivaara & Väinö-Waltteri Granat

# Signal processing in digital holography - Report

# **Abstract**

Roosa Kuusivaara & Väinö-Waltteri Granat: Signal processing in digital holography  
- Report

Laboratory Report

Tampere University

Master's Degree Programme in Signal Processing

November 2023

---

This report documents the work done in the Signal processing in digital holography assignment as a part of the Advanced signal processing laboratory course. In this assignment we familiarized ourselves with basics of interference based holography.

**Keywords:** Hologram, Mach-Zehnder interferometer, Wave propagation.

# Contents

1	Introduction . . . . .	1
2	Methodology . . . . .	2
2.1	Experiment setup . . . . .	2
2.2	Object reconstruction . . . . .	2
3	Results . . . . .	5
3.1	Wavefront Reconstruction . . . . .	5
3.2	Analysis . . . . .	5
4	Conclusions . . . . .	7
	References . . . . .	8
	APPENDIX . . . . .	8

## 1 Introduction

In this report we describe our work done in the 'Signal processing in digital holography' laboratory assignment for the Advanced Signal Processing Laboratory course. In this assignment we used a Mach-Zehnder interferometer to record holographic images and Matlab to reconstruct the object from the recorded wavefront.

## 2 Methodology

This section describes the method we used to capture the holograms and reconstruct the objects.

### 2.1 Experiment setup

Holographic imagery relies on the phenomena of light wave superpositions. In this assignment the object is captured using a Mach-Zehnder interferometer, where a monochromatic laser is used to produce two wavefronts. The first wavefront passes through the object which we want to record, and the second wavefront is used as reference wavefront. The wavefronts are combined to a single wavefront via the interference phenomena, so that one of the wavefronts is combined with the other wavefront in a slight angle, causing a phase difference. The resulting total wave is then captured with a CMOS camera to produce an image file. The phase difference is detectable as fringes in the captured images. The complete setup for the experiment is shown in figure 2.1.


### 2.2 Object reconstruction

The object can be extracted from the captured hologram, using a two step method.


First step is to use Fourier filtering to extract the object wave. Since we know that the hologram is defined by the following equation:

$$H(x, y) = E_0(x, y)^2 + E_r(x, y)^2 + U_0(x, y)U_r^*(x, y) + U_0^*(x, y)U_r(x, y) \quad (2.1)$$

Where  $H(x, y)$  is the hologram wavefront in given position,  $E_0$  is the objects

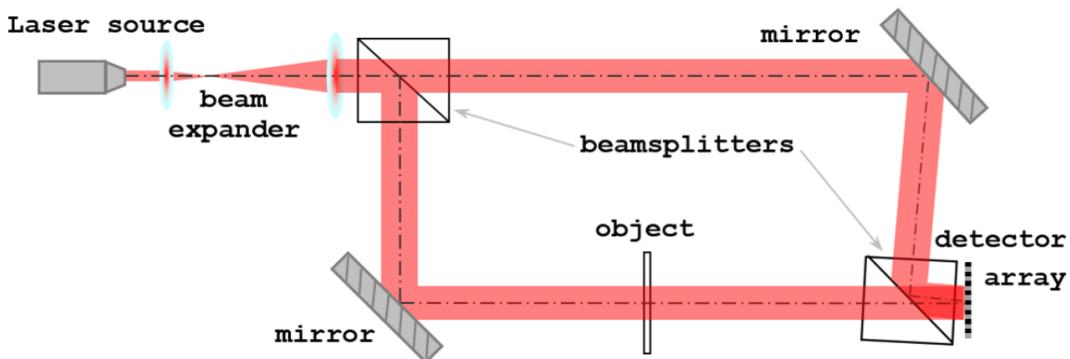

**Figure 2.1** Experiment was done using a Mach-Zehnder interferometer. Taken from [assignment]

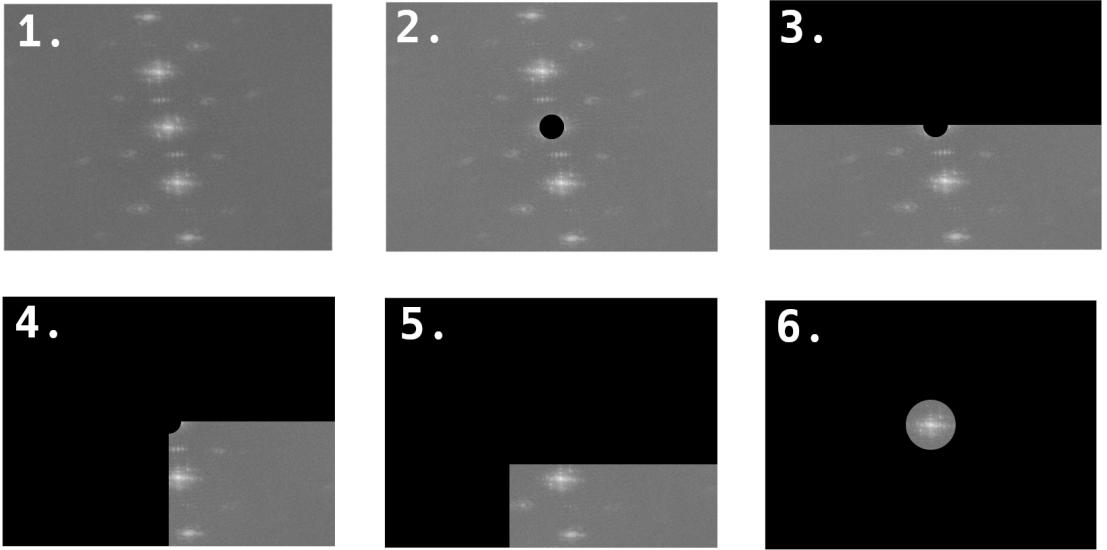

**Figure 2.2** Extracting the second term from Fourier transformation of the hologram

amplitude,  $E_r$  is the reference wave's amplitude,  $U_r$  is the reference wavefront and  $U_0$  is the object wavefront which we want to extract. The extraction can be accomplished by performing frequency filtering to the Fourier transformation of the hologram, which is shown in equation 2.2. From the Fourier transformation we can see that the second term of the hologram contains the object wavefront  $U_0$ . Equations 2.3 and 2.4 show the Fourier transformations for the second and third term respectively.

$$\begin{aligned} \mathcal{F}[H(x, y)] = & \mathcal{F}[E_0(x, y)^2 + E_r(x, y)^2] + \mathcal{F}[U_0(x, y)E_r \exp(i2\pi\eta)] \\ & + \mathcal{F}[U_0^*(x, y)E_r \exp(-i2\pi\eta)] \end{aligned} \quad (2.2)$$

$$F[U_0(x, y)E_r \exp(i2\pi\eta)] = F[U_0]E_r \otimes \delta(f - \eta) \quad (2.3)$$

$$F[U_0^*(x, y)E_r \exp(-i2\pi\eta)] = F[U_0^*]E_r \otimes \delta(f + \eta) \quad (2.4)$$

By removing all the other terms from the Fourier transformation of the hologram, we are left with only the Fourier transformation of the object wave. The object wave can be brought back to spatial domain with inverse Fourier transformation. The figure 2.2 shows how the second term is isolated from the rest of the hologram.

The wavefront alone isn't that usable. The object needs to be projected to an object plane using wave propagation. Using Rayleigh-Sommerfield holographic model with the angular spectrum transfer function defined in equation 2.5 the

object is projected on the object plane as defined by the equation 2.6, where  $z$  is the propagation distance. To reconstruct an object that is in focus, we need to determine the optimal distance for the wave to be propagated. This can be accomplished with an autofocusing algorithm or manually by comparing reconstruction results with different values of  $z$  by eye.

$$H(f_x, f_y, z) = \begin{cases} \exp[i2\pi z \sqrt{1 - \lambda^2(f_x^2 + f_y^2)}], & f_x^2 + f_y^2 \leq 1/\lambda^2 \\ 0, & \text{otherwise} \end{cases} \quad (2.5)$$

$$U_0(x, y, z) = F^{-1}[H(f_x, f_y, z) \cdot F[U_0(x, y, 0)]] \quad (2.6)$$

## 3 Results


This section of report presents the results of the wavefront reconstruction, and the analysis of those results. The object wavefront reconstruction process described in the section 2 was implemented as a Matlab script, using a provided implementation of the wave propagation algorithm.

### 3.1 Wavefront Reconstruction


In the laboratory we captured two hologram containing the same object with two different fringe widths. The distance between ray and camera was 11 cm, wavelength was  $532 * 10^{-9}$  m, and a number of pixels was  $3,45 * 10^{-6}$  m. The first hologram had fringes of 5 pixels and the second had fringes of 12 pixels.

### 3.2 Analysis

The figures 3.1 and 3.2 show the amplitude and phase images of the reconstructed object wavefronts. In 5 pixel fringe hologram, for getting the sharpest details of phase, the best propagation distance was 29,77 cm. In 12 pixel fringe hologram, the sharpest details of phase was determined to be 29,86 cm.

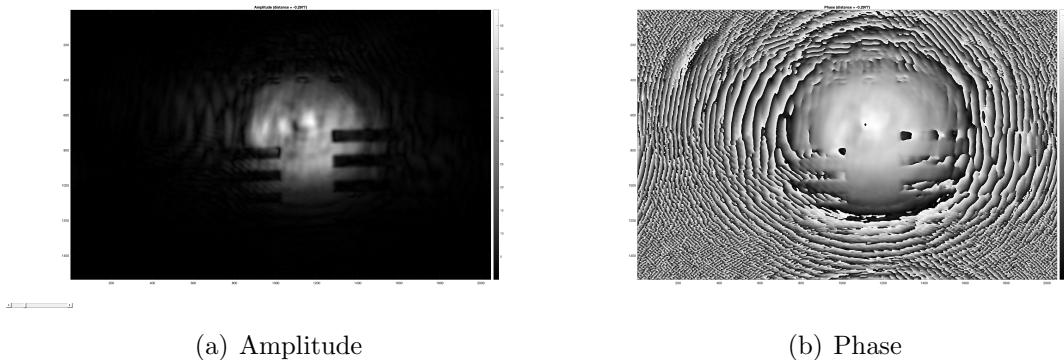


**Figure 3.1** Object wavefront reconstruction of the 5px fringe image

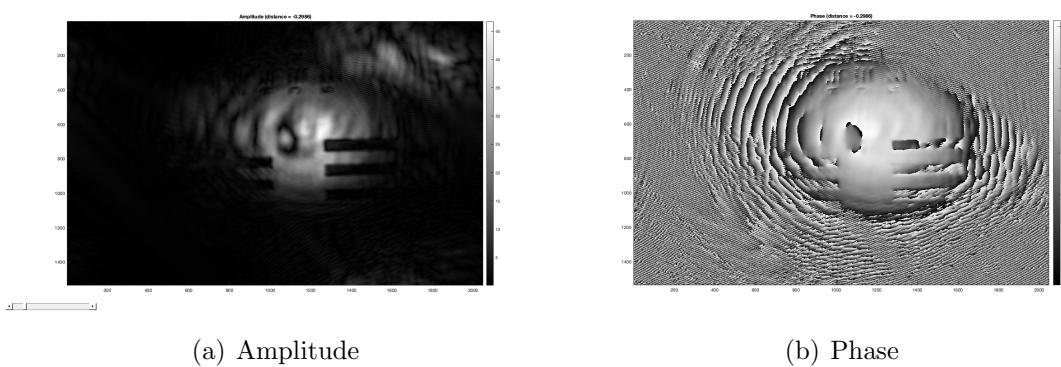


**Figure 3.2** Object wavefront reconstruction of the 12px fringe image

## 4 Conclusions

## APPENDIX

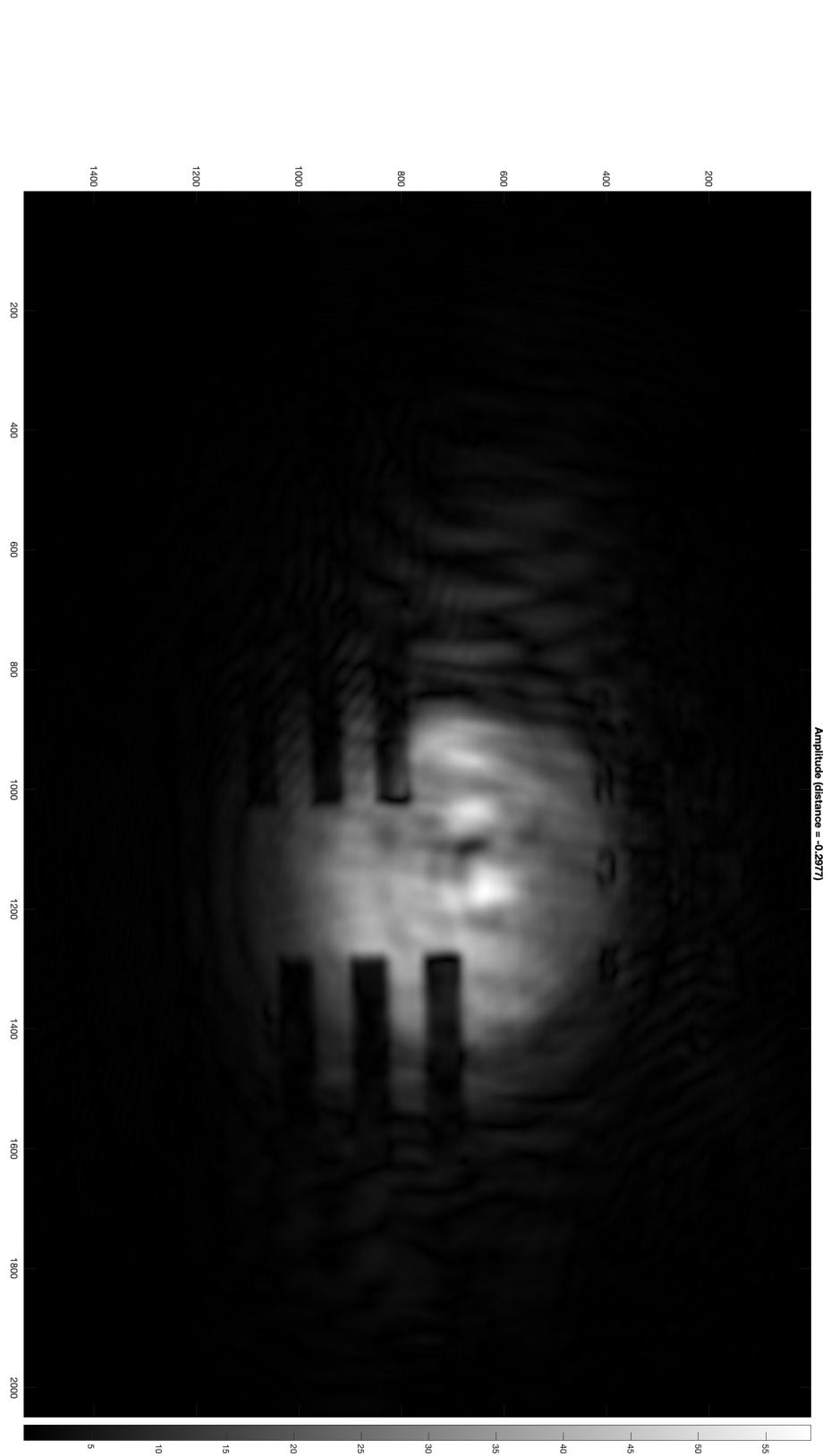

**Figure 1** Caption for your image

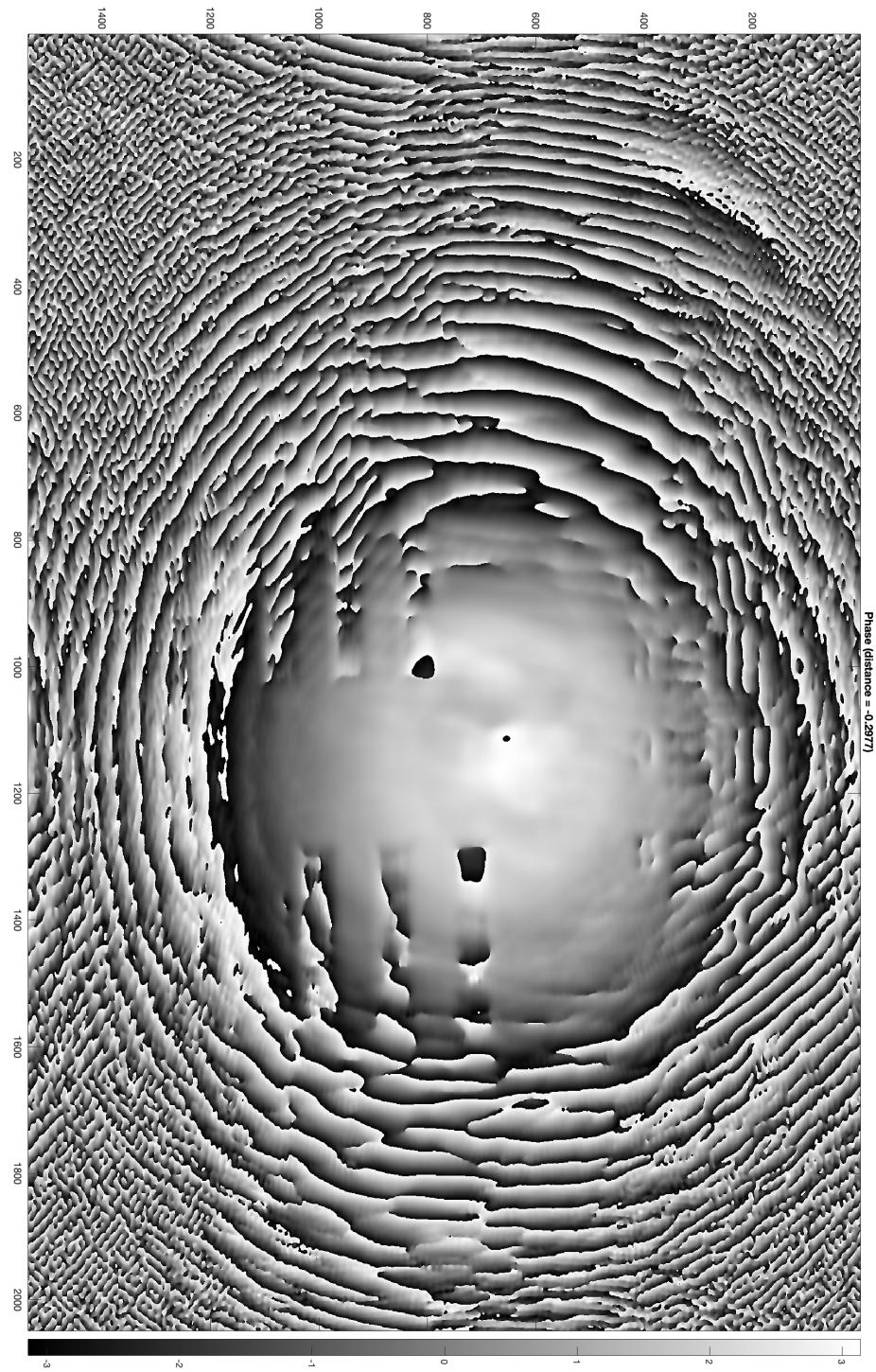

**Figure 2** Caption for your image